

Adult ED & FSED IV Vancomycin Update- Go Live 12/13/2022




To better improve administration of the full dose, the medication order for **Vancomycin Doses greater than 1500mg** was updated. The medication will be profiled simultaneously and will display as two entries within eMAR, scheduled a few minutes apart. This will allow the Nurse to scan each entry and to program the pump to **infuse concurrently** as per administration instructions found within the label comments.

Vancomycin infusing via Concurrent Mode:

NOTE: Utilize Smart Pump Programming as applicable for your facility

To Manually Program:

- Turn on IV Pump and select the appropriate **CCA (ER Adults)**
- In eMAR:
 - Scan the **patient**
 - Scan **bag #1**
 - Select the **appropriate medication** in the Drug Library
 - Program bag on **channel A per Label Comments**
 - Press the **Start** button and **confirm** to start infusion
 - Scan **bag #2**
 - Select the **appropriate medication** in the Drug Library
 - Program bag on **channel B per Label Comments**
 - When programming on channel B, the *Mode* will display as “Piggyback”
 - Press the **Change Mode option**, it will now display “Concurrent”
 - Press the **Start** button and **confirm** to start infusion

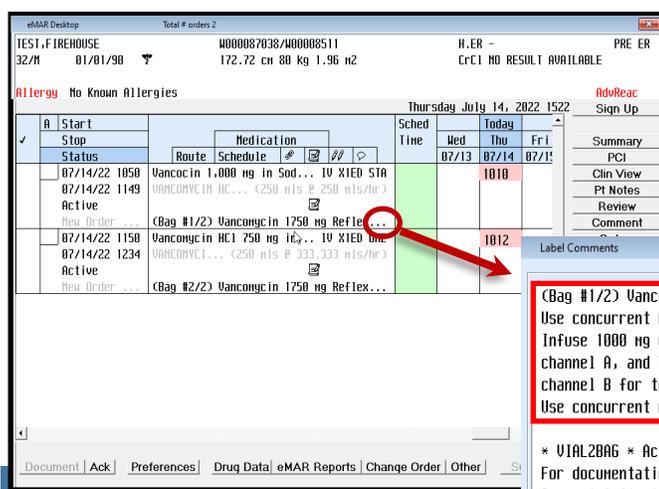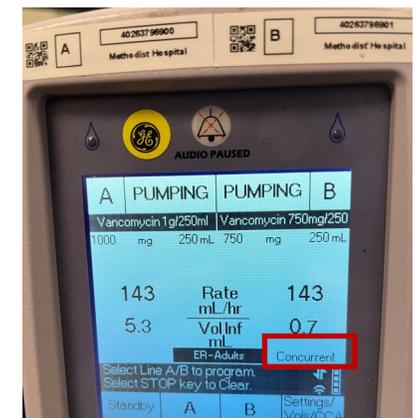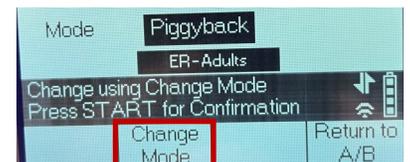

How to view Label Comments

Click on the Ellipsis (...)

To display label comments, which contain important instructions to follow for medication administration

